



Infrared-Thermometer 31.1133.01
for temperature measurement on
the forehead




GB

Quick Guide

1. Download the app OUcare (IRT OUcare) from the Apple App Store for your iPhone or from the Google Play store for your Android smartphone or open the already downloaded app.



Launch the app IRT OUcare

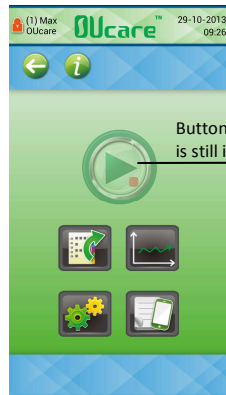

Button
is still inactive

2. Place the O/S switch in the correct position of the operating system of your smartphone (switch position 1 for Android, 2 Apple iOS for).



Protective cap
Schutzkappe

3. Plug the device into the earphone jack of your smartphone and remove the protective cap.




Button becomes active

- Turn it to a suitable position for your measurement.



4. Direct the sensor to the forehead. Keep a distance of 1 to 3 cm to the forehead! The measuring time is approximately 1 second. Make sure that the distance to the forehead is not changed during the measurement.

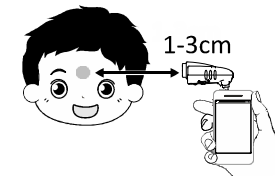

5. The last user (1-9) will appear in the top left corner. To start the measurement press the round button on the screen.




With activated voice output you will hear the measured temperature. The unit is now ready for another test. (Display L/H: temperature is out of range)

Temperature measurement of an object or surface

By touching the button



you change the mode to the measure of the temperature of an object or surface.

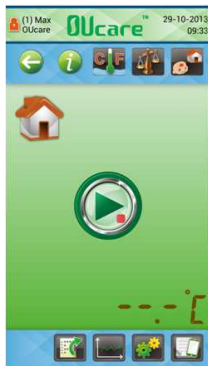

The procedure is the same as for the measurement on the forehead. But there is no voice output.

Additional functions of the IRT OUCare app

Change between °C / °F

- By touching the button



the values are displayed in the selected unit.

Offset-Kompensation

- By touching the button



an adjustment to any future measurement can be set. This function is only intended for the measurement of body temperature, e.g. to adjust the individual normal temperature to 37 °C.

User defined settings

- By touching the button



you will enter the language setting.

Specific languages, available in the paid version only.



Voice output on/off



- By touching the button



you will be able to set user profiles (max. 9 user). Select a specific user.

- By touching the button



you get to set up to 5 email addresses to send your measurements (1st line: name, 2nd line: email address).

- By touching the button



you are able to set up to 5 mobile numbers to send your measurements

via SMS service- (1st line: name, 2nd line: mobile number).

- By touching the button



the current settings are displayed. Measuring object, temperature unit and offset compensation (at body-temperature) can be changed.

Graphic display of measured values

- By touching the button



you see the measured values in a chart. T = individual values, D = daily average. Changing between object and body temperature is possible.

List of measured recorded values

- By touching the button



you get a list of history readings in chronological order. Change between


object and body temperature, display, sending and deleting of single values possible.

- By touching the button



you are able to add a value manually.

The LED-Indicator

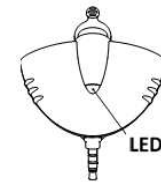

red: no contact or empty battery

green: Function OK

CAUTION

Disconnect the infrared thermometer after use from your smartphone. Only then will it turn off automatically. Otherwise, it will still consume power.